Technical Data Sheet

FITC Mouse Anti-Human CD61

Product Information

Material Number: 561913

Alternate Name: Integrin β3 Chain

25 tests Size Vol. per Test: 20 µl VI-PL2 Clone: Mouse IgG1, κ Isotype: Reactivity: QC Testing: Human


Tested in Development: Dog

Workshop: V 5T-124

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Reacts with integrin $\beta 3$, a 105 kDa transmembrane glycoprotein, expressed on platelets, megakaryocytes, osteoclasts and endothelia. Integrin β3 associates with gpIIa (CD41) to form the CD41/CD61 complex which mediates platelet adhesion and aggregation. CD61 also associates with CD51 to form the CD51/CD61 complex (vitronectin receptor). CD61 appears to bind to fibrinogen, fibronectin, vWF, vitronectin, and thrombospondin to mediate cell adhesion.

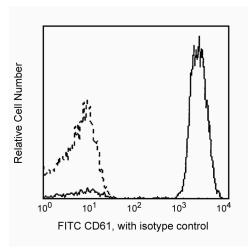

Profile of peripheral blood platelets analyzed on a BD FACScan™ (BDIS, San Jose, CA)

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone
555748	FITC Mouse IgG1, κ Isotype Control	100 tests	MOPC-21
554656	Stain Buffer (FBS)	500 ml	(none)

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use $1 \times 10^{\circ}6$ cells in a 100- μ l experimental
- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

BD Biosciences

bdbiosciences.com

United States Asia Pacific Latin America/Caribbean 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



561913 Rev. 1

- 4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 5. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 6. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 7. An isotype control should be used at the same concentration as the antibody of interest.

References

Schlossman S, Boumell L, et al, ed. Leucocyte Typing V. New York: Oxford University Press; 1995. (Biology)

561913 Rev. 1 Page 2 of 2