

Technical Data Sheet


BV421 Mouse Anti-Human CD147**Product Information**

| | |
|-------------------------|---|
| Material Number: | 562583 |
| Alternate Name: | BSG; Basigin; Basi; Neurothelin; 5F7; EMMPRIN; M6; OK; TCSF |
| Size: | 50 tests |
| Vol. per Test: | 5 µl |
| Clone: | HIM6 |
| Isotype: | Mouse IgG1, κ |
| Reactivity: | QC Testing: Human |
| Workshop: | VI NL109 |
| Storage Buffer: | Aqueous buffered solution containing BSA and ≤0.09% sodium azide. |

Description

The HIM6 monoclonal antibody specifically binds to CD147 which is encoded by BSG. CD147 is a type I transmembrane glycoprotein (30-50 kDa) of the immunoglobulin super-gene family. Neurothelin, a blood-brain barrier-specific molecule, was clustered as CD147 in the Sixth Human Leukocyte Differentiation Antigen (HLDA) workshop. It bears homology with mouse gp42 or basigin, human "M6" or "EMMPRIN", rat OX-47 or CD-9, and avian HT7 or 5A11. CD147 is also known as Tumor cell-derived collagenase stimulatory factor (TCSF). CD147 is a molecule that is broadly expressed on cells of hematopoietic and non-hematopoietic origin. Its expression on specific cell types may be regulated by cytokines. CD147 plays a role in embryonal blood-brain barrier development and a role in integrin-mediated adhesion in brain endothelia.

The antibody was conjugated to BD Horizon™ BV421 which is part of the BD Horizon™ Brilliant Violet™ family of dyes. With an Ex Max of 407-nm and Em Max at 421-nm, BD Horizon™ BV421 can be excited by the violet laser and detected in the standard Pacific Blue™ filter set (eg, 450/50-nm filter). BD Horizon™ BV421 conjugates are very bright, often exhibiting a 10 fold improvement in brightness compared to Pacific Blue™ conjugates.

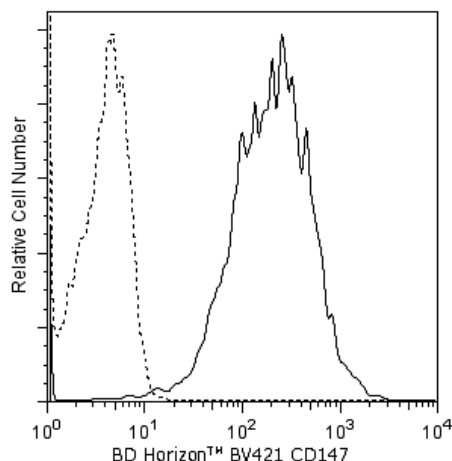

Flow cytometric analysis of CD147 expression on human peripheral blood lymphocytes. Whole blood was stained with BD Horizon™ BV421 Mouse Anti-Human CD147 antibody (Cat. No. 562583; solid line histogram), or with a BD Horizon™ BV421 Mouse IgG1, κ Isotype Control (Cat. No. 562438; dashed line histogram). Erythrocytes were lysed with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from gated events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV421 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV421 were removed.

Application Notes**Application**

| | |
|----------------|------------------|
| Flow cytometry | Routinely Tested |
|----------------|------------------|

BD Biosciences

bdbiosciences.com

| United States | Canada | Europe | Japan | Asia Pacific | Latin America/Caribbean |
|---------------|--------------|---------------|--------------|--------------|-------------------------|
| 877.232.8995 | 800.979.9408 | 32.53.720.550 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995 |

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Suggested Companion Products

| Catalog Number | Name | Size | Clone |
|----------------|-------------------------------------|--------|--------|
| 562438 | BV421 Mouse IgG1, k Isotype Control | 50 µg | X40 |
| 555899 | Lysing Buffer | 100 ml | (none) |
| 554656 | Stain Buffer (FBS) | 500 ml | (none) |

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
3. An isotype control should be used at the same concentration as the antibody of interest.
4. Please refer to www.bdbiosciences.com/pharming/en/protocols for technical protocols.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. Pacific Blue™ is a trademark of Molecular Probes, Inc., Eugene, OR.
8. Brilliant Violet™ 421 is a trademark of Sirigen.

References

Majdic O, Pickl WF, Kohl P, Stockinger H, Knapp W. EC16.3 CD147 Workshop: Reactivity and epitope mapping of CD147 monoclonal antibodies. In: Kishimoto T, Kikutani H, von dem Borne AEGK, ed. *Leukocyte Typing VI: White Cell Differentiation Antigens*. New York: Garland Publishing Inc; 1998:765-766. (Clone-specific: Flow cytometry)

Riethdorf S, Reimers N, Assmann V, Kornfeld JW, Terracciano L, Sauter G, Pantel K. High incidence of EMMPRIN expression in human tumors. *Int J Cancer*. 2006; 119(8):1800-1810. (Clone-specific: Immunohistochemistry, Immunoprecipitation, Western blot)

Rizzo A, Aragona E, Dino O, et al. EC16.1 CD147 Workshop: Expression of CD147 (neurothelin) in liver and lung cancer. In: Kishimoto T, Kikutani H, von dem Borne AEGK, ed. *Leukocyte Typing VI: White Cell Differentiation Antigens*. New York: Garland Publishing Inc; 1998:763-764. (Clone-specific: Immunohistochemistry)

Stockinger H, Ebel T, Hansmann C, et al. EC16 CD147 (neurothelin/basigin) Workshop Panel Report. In: Kishimoto T, Kikutani H, von dem Borne AEGK, ed. *Leukocyte Typing VI: White Cell Differentiation Antigens*. New York: Garland Publishing Inc; 1998:760-763. (Clone-specific: Flow cytometry, Immunoprecipitation)

Sudou A, Ozawa M, Muramatsu T. Lewis X structure increases cell substratum adhesion in L cells. *J Biochem (Tokyo)*. 1995; 117(2):271-275. (Biology)

BD Biosciences

bdbiosciences.com

| | | | | | |
|---------------|--------------|---------------|--------------|--------------|-------------------------|
| United States | Canada | Europe | Japan | Asia Pacific | Latin America/Caribbean |
| 877.232.8995 | 800.979.9408 | 32.53.720.550 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995 |

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

