# Technical Data Sheet

# BV605 Mouse Anti-Human CD117

### **Product Information**

**Material Number:** 562687

Alternate Name: KIT; c-Kit; SCFR; PBT; Mast/stem cell growth factor receptor

50 Tests Size Vol. per Test: 5 μl 104D2 Clone:

Megakaryoctic cell line MOLM-1 Immunogen:

Isotype: Mouse (BALB/c) IgG1 Reactivity: QC Testing: Human

VI C30 Workshop:


Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

## Description

The 104D2 monoclonal antibody specifically binds to human CD117, the receptor for stem cell factor (SCF). It selectively recognizes NIH-3T3 cells transfected with human c-kit, the gene that codes for SCF-R. The 104D2 antibody does not block the epitope that binds SCF. In the bone marrow of humans and mice, SCF is expressed primarily on hematopoietic progenitor cells. Lack of functional SCF or deficient SCF-R caused by mutations in the SI and W loci, respectively, can result in severe anemia and a decrease in the number of primitive progenitor cells in mice. Human hematopoietic progenitor cells can be recognized by their surface expression of CD34. This cell population constitutes a small subset (1% to 5%) of bone marrow cells. CD34+ cells contain a small subpopulation of primitive/non-committed progenitors, with the remaining fraction being cells committed to the various hematopoietic lineages. SCF alone induces extensive proliferation of erythroid-committed progenitor cells (CD34lo CD71hi CD64-). On primitive (CD34hi CD38lo CD50+) and granulo-monocytic (CD34+ CD64+) progenitor cells, SCF synergistically enhances the effects of other cytokines, the strongest of which are on the primitive progenitor cells. In addition, SCF promotes survival of primitive progenitors in the absence of proliferation. The receptor is highly expressed at similar levels on all of the three mentioned CD34+ cell subsets, whereas B-lymphoid committed progenitor cells (CD34+ CD19+) express low levels of SCF-R. Among CD34- bone marrow cells, only a small number of cells (mostly erythroid) express the receptor.

This antibody is conjugated to BD Horizon BV605 which is part of the BD Horizon Brilliant™ Violet family of dyes. With an Ex Max of 407-nm and Em Max of 602-nm, BD Horizon BV605 can be excited by a violet laser and detected with a standard 610/20-nm filter set. BD Horizon BV605 is a tandem fluorochrome of BD Horizon BV421 and an acceptor dye with an Em max at 605-nm. Due to the excitation of the acceptor dye by the green (532 nm) and yellow-green (561 nm) lasers, there will be significant spillover into the PE and BD Horizon PE-CF594 detectors off the green or yellow-green lasers. BD Horizon BV605 conjugates are very bright, often exhibiting brightness equivalent to PE conjugates and can be used as a third color off of the violet laser.

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794).

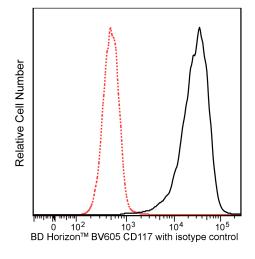

Flow cytometric analysis of CD117 expression on human TF-1 Cells. TF-1 cells (Human ervthroleukemia cell line; ATCC Cat. No. CRL-2003) were stained with the BD Horizon™ BV605 Mouse Anti-Human CD117 antibody (Cat. No. 562687; solid line histogram) or with BD Horizon<sup>T</sup> BV605 Mouse IgG1, κ Isotype Control (Cat. No. 562652; dashed line histogram). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of viable cells. Flow cytometry was performed using a BD $^{\text{TM}}$  LSR II Flow Cytometry System.

## **BD Biosciences**

bdbiosciences.com

 
 Canada
 Europe
 Japan

 800.268.5430
 32.2.400.98.95
 0120.8555.90
 **United States** Asia Pacific Latin America/Caribbean

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale. Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



### Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV605 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV605 were removed.

### **Application Notes**

## Application

Flow cytometry	Routinely Tested	
1 low cytometry	Routinery rested	

## **Suggested Companion Products**

Catalog Number	Name	Size	Clone
562652	BV605 Mouse IgG1, κ Isotype Control	50 μg	X40
554656	Stain Buffer (FBS)	500 mL	(none)
563794	Brilliant Stain Buffer	5 mL	(none)

### **Product Notices**

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10<sup>6</sup> cells in a 100-μl experimental
- An isotype control should be used at the same concentration as the antibody of interest.
- Although every effort is made to minimize the lot-to-lot variation in the efficiency of the fluorochrome energy transfer, differences in the residual emission from BD Horizon<sup>TM</sup> BV421 may be observed. Therefore, we recommend that individual compensation controls be performed for every BD Horizon<sup>TM</sup> BV605 conjugate.
- Please observe the following precautions: Absorption of visible light can significantly alter the energy transfer occurring in any tandem fluorochrome conjugate; therefore, we recommend that special precautions be taken (such as wrapping vials, tubes, or racks in aluminum foil) to prevent exposure of conjugated reagents, including cells stained with those reagents, to room illumination.
- Texas Red is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States. 6.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 10. CFTM is a trademark of Biotium, Inc.

#### References

Ashman LK, Buhring HJ, Aylett GW, Broudy VC, Muller C. Epitope mapping and functional studies with three monoclonal antibodies to the c-kit receptor tyrosine kinase, YB5.B8, 17F11, and SR-1. J Cell Physiol. 1994; 158(3):545-554. (Clone-specific: Flow cytometry)

Berardi AC, Wang A, Levine JD, Lopez P, Scadden DT. Functional isolation and characterization of human hematopoietic stem cells. Science. 1995; 267:104-108. (Clone-specific)

Bernstein ID, Andrews RG, Zsebo KM. Recombinant human stem cell factor enhances the formation of colonies by CD34+ and CD34+lin- cells, and the generation of colony-forming cell progeny from CD34+lin- cells cultured with interleukin-3, granulocyte colony-stimulating factor, or granulocyte-macrophage colony-stimulating factor. Blood. 1991; 77:2316-2321. (Biology)

Briddell RA, Broudy VC, Bruno E, Brandt JE, Srour EF, Hoffman R. Further phenotypic characterization and isolation of human hematopoietic progenitor cells using a monoclonal antibody to the c-kit receptor. Blood. 1992; 79:3159-3167. (Clone-specific)

Broudy VC, Lin N, Zsebo KM, Birkett NC, Smith KA, Bernstein ID, Papayannopoulou T. Isolation and characterization of a monoclonal antibody that recognizes the human c-kit receptor, Blood, 1992; 79:338-346, (Clone-specific)

de Vries P, Brasel KA, Eisenman JR, Alpert AR, Williams DE. The effect of recombinant mast cell growth factor on purified murine hematopoietic stem cells. J Exp Med. 1991; 173:1205-1211. (Biology)

Gunji Y, Nakamura M, Osawa H, Naqayoshi K, Nakauchi H, Miura Y, Yanagisawa M, Suda T. Human primitive hematopoietic progenitor cells are more enriched in KITlow cells than in KIThigh cells. Blood. 1993; 82:3283-3289. (Clone-specific) Ikuta K, Weissman IL. Evidence that hematopoietic stem cells express mouse c-kit but do not depend on steel factor for their generation. Proc Natl Acad Sci U S

A. 1992; 89(4):1502-1506. (Clone-specific) Keller JR, Ortiz M, Ruscetti FW. Steel factor (c-kit ligand) promotes the survival of hematopoietic stem/progenitor cells in the absence of cell division. Blood. 1995;

86:1757-1764. (Biology)

Kinashi T, Springer TA. Steel factor and c-kit regulate cell-matrix adhesion. Blood. 1994; 83:1033-1038. (Biology)

Lansdorp PM, Dragowska W, Long-term erythropoiesis from constant numbers of CD34+ cells in serum-free cultures initiated with highly purified progenitor cells from human bone marrow. J Exp Med. 1992; 175:1501-1509. (Biology)

Matsuo Y, Adachi T, Tsubota T, Imanishi J, Minowada J. Establishment and characterization of a novel megakaryoblastic cell line, MOLM-1, from a patient with chronic myelogenous leukemia. Hum Immunol. 1991; 4:261-264. (Biology)

Okada S, Nakauchi H, Nagayoshi K, Nishikawa S, Miura Y, Suda T. Enrichment and characterization of murine hematopoietic stem cells that express c-kit molecule. Blood. 1991; 78:1706-1712. (Clone-specific)

Olweus J, Lund-Johansen F, Terstappen LW. CD64/Fc gamma RI is a granulo-monocytic lineage marker on CD34+ hematopoietic progenitor cells. Blood. 1995;

Olweus J, Terstappen LW, Thompson PA, Lund-Johansen F. Expression and function of receptors for stem cell factor and erythropoietin during lineage commitment of human hematopoietic progenitor cells. Blood. 1996; 88:1594-1607. (Clone-specific)

Ratajczak MZ, Luger SM, DeRiel K, Abrahm J, Calabretta B, Gewirtz AM. Role of the KIT protooncogene in normal and malignant human hematopoiesis. Proc Natl Acad Sci U S A. 1992; 89:1710-1714. (Biology)

Simmons PJ, Aylett GW, Niutta S, To LB, Juttner CA, Ashman LK. c-kit is expressed by primitive human hematopoietic cells that give rise to colony-forming cells in stroma-dependent or cytokine-supplemented culture. Exp Hematol. 1994; 22:157-165. (Clone-specific)

## **BD Biosciences**

bdbiosciences.com

 
 Canada
 Europe
 Japan

 800.268.5430
 32.2.400.98.95
 0120.8555.90
 **United States** Asia Pacific Latin America/Caribbean

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violatio of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with thuse of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.
For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



Page 2 of 3

Terstappen LW, Huang S, Safford M, Lansdorp PM, Loken MR. Sequential generations of hematopoietic colonies derived from single nonlineage-committed CD34+CD38- progenitor cells. *Blood.* 1991; 77(6):1218-1227. (Biology)

Waller EK, Olweus J, Lund-Johansen F, Huang S, Nguyen M, Guo GR, Terstappen L. The "common stem cell" hypothesis reevaluated: human fetal bone marrow contains separate populations of hematopoietic and stromal progenitors. Blood. 1995; 85:2422-2435. (Biology)

Yarden Y, Kuang WJ, Yang-Feng T, et al. Human proto-oncogene c-kit: a new cell surface receptor tyrosine kinase for an unidentified ligand. EMBO J. 1987; 6(11):3341-3351. (Clone-specific)

# **BD Biosciences**

bdbiosciences.com

 
 United States
 Canada
 Europe
 Japan

 877.232.8995
 800.268.5430
 32.2.400.98.95
 0120.8555.90
 **Asia Pacific** 65.6861.0633 Latin America/Caribbean

For country contact information, visit **bdbiosciences.com/contact** 

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



562687 Rev. 3 Page 3 of 3