# **Technical Data Sheet**

# PE-CF594 Rat Anti-Mouse Ly-6G

#### **Product Information**

Material Number: 562700

Alternate Name: Ly6g; Lymphocyte antigen 6G; Lymphocyte antigen 6 complex, locus G; Gr1

 Size:
 50 μg

 Concentration:
 0.2 mg/ml

 Clone:
 1A8

Immunogen: Ly-6G-transfected EL4J cell line


Isotype:Rat (LEW) IgG2a,  $\kappa$ Reactivity:QC Testing: Mouse


Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

#### Description

The 1A8 antibody reacts with Ly-6G, a 21-25-kDa GPI-anchored protein. In the bone marrow, Ly-6G is expressed on the majority of the largest cells, which are predominantly granulocytes, and not on lymphoid or erythroid cells. In the periphery, it is expressed on granulocytes. The mAb RB6-8C5 (Cat. No. 557445/553123) recognizes both Ly-6G and Ly-6C and blocks the binding of mAb 1A8 to Ly-6G.

This antibody is conjugated to BD Horizon<sup>TM</sup> PE-CF594, which has been developed exclusively by BD Biosciences as a better alternative to PE-Texas Red®. PE-CF594 excites and emits at similar wavelengths to PE-Texas Red® yet exhibits improved brightness and spectral characteristics. Due to PE having maximal absorption peaks at 496 nm and 564 nm, PE-CF594 can be excited by the blue (488-nm), green (532-nm) and yellow-green (561-nm) lasers and can be detected with the same filter set as PE-Texas Red® (eg 610/20-nm filter).

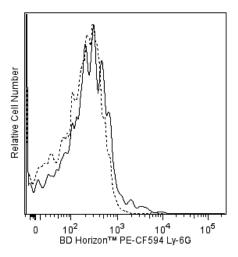

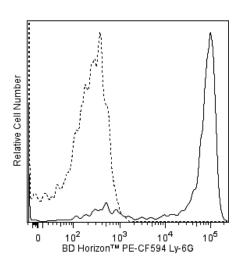

Flow cytometric analysis of Ly-6G expression on mouse bone-marrow leukocytes. A BALB/c bone-marrow cell suspension was stained with either BD Horizon™ PE-CF594 Rat Anti-Mouse Ly-6G antibody (Cat. No. 562700; solid line histogram) or BD Horizon™ PE-CF594 Rat IgG2a, κ Isotype Control (Cat. No. 562302; dashed line histogram). Flow cytometric histograms were derived from gated events based on the light scattering characteristics of viable lymphoid (Left Panel) or myeloid (Right Panel) cells. Flow cytometry was performed using a BD™ LSR II Flow Cytometry System.

### **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The antibody was conjugated with BD Horizon™ PE-CF594 under optimum conditions, and unconjugated antibody and free PE-CF594 were removed.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

#### **Application Notes**

Application

Flow cytometry Routinely Tested

#### **BD Biosciences**

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 800.979.9408
 32.53.720.550
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



#### **Suggested Companion Products**

| Catalog Number | Name                                  | Size   | Clone  |
|----------------|---------------------------------------|--------|--------|
| 562302         | PE-CF594 Rat IgG2a, κ Isotype Control | 0.1 mg | R35-95 |
| 554656         | Stain Buffer (FBS)                    | 500 ml | (none) |

#### **Product Notices**

- 1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 3. An isotype control should be used at the same concentration as the antibody of interest.
- 4. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 5. Please observe the following precautions: Absorption of visible light can significantly alter the energy transfer occurring in any tandem fluorochrome conjugate; therefore, we recommend that special precautions be taken (such as wrapping vials, tubes, or racks in aluminum foil) to prevent exposure of conjugated reagents, including cells stained with those reagents, to room illumination.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 7. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- 8. Texas Red is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- 9. CF<sup>TM</sup> is a trademark of Biotium, Inc.
- 10. When excited by the yellow-green (561-nm) laser, the fluorescence may be brighter than when excited by the blue (488-nm) laser.
- 11. This product is provided under an Agreement between BIOTIUM and BD Biosciences. The manufacture, use, sale, offer for sale, or import of this product is subject to one or more patents or pending applications owned or licensed by Biotium, Inc. This product, and only in the amount purchased by buyer, may be used solely for buyer's own internal research, in a manner consistent with the accompanying product literature. No other right to use, sell or otherwise transfer (a) this product, or (b) its components is hereby granted expressly, by implication or by estoppel. This product is for research use only. Diagnostic uses require a separate license from Biotium, Inc. For information on purchasing a license to this product including for purposes other than research, contact Biotium, Inc., 3159 Corporate Place, Hayward, CA 94545, Tel: (510) 265-1027. Fax: (510) 265-1352. Email: btinfo@biotium.com.
- Because of the broad absorption spectrum of the tandem fluorochrome, extra care must be taken when using multi-laser cytometers, which
  may directly excite both PE and CFTM594.

#### References

Fleming TJ, Fleming ML, Malek TR. Selective expression of Ly-6G on myeloid lineage cells in mouse bone marrow. RB6-8C5 mAb to granulocyte-differentiation antigen (Gr-1) detects members of the Ly-6 family. *J Immunol.* 1993; 151(5):2399-2408. (Immunogen: Flow cytometry, Immunoprecipitation)
Fleming TJ, Malek TR. Multiple glycosylphosphatidylinositol-anchored Ly-6 molecules and transmembrane Ly-6E mediate inhibition of IL-2 production. *J Immunol.* 1994; 153(5):1955-1962. (Biology)

## **BD Biosciences**

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 800.979.9408
 32.53.720.550
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale. Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562700 Rev. 1 Page 2 of 2