

## Technical Data Sheet


**BV510 Mouse Anti-Human CD45RA****Product Information**

<b>Material Number:</b>	<b>563031</b>
<b>Alternate Name:</b>	CD45R; PTPRC; LCA; Leukocyte common antigen
<b>Size:</b>	50 tests
<b>Vol. per Test:</b>	5 µl
<b>Clone:</b>	HI100
<b>Isotype:</b>	Mouse IgG2b, κ
<b>Reactivity:</b>	QC Testing: Human
<b>Workshop:</b>	IV N906
<b>Storage Buffer:</b>	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

**Description**

The HI100 monoclonal antibody specifically binds to the 220 kDa isoform of the human leukocyte common antigen, CD45RA. CD45RA is expressed on approximately 40-50% of peripheral CD4+ T cells, 50% of peripheral CD8+ T cells and on a portion of B cells and monocytes. The CD45RA antigen is expressed by naïve and activated T cells. CD45RA-specific antibodies are useful for the study of the suppressor/inducer subpopulation of CD4+ lymphocytes.

The antibody was conjugated to BD Horizon™ BV510 which is part of the BD Horizon™ Brilliant Violet™ family of dyes. With an Ex Max of 405-nm and Em Max at 510-nm, BD Horizon™ BV510 can be excited by the violet laser and detected in the BD Horizon V500 (525/50-nm) filter set. BD Horizon™ BV510 conjugates are useful for the detection of dim markers off the violet laser.

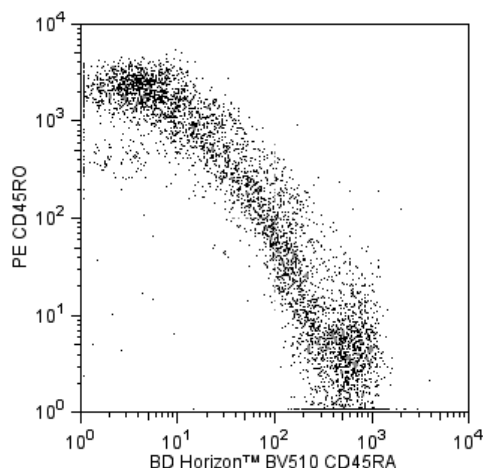

**Two-color flow cytometric analysis of CD45RA expression on human peripheral blood lymphocytes.** Human whole blood was stained with BD Horizon™ BV510 Mouse Anti-Human CD45RA (Cat. No. 563031) and PE Mouse Anti-Human CD45RO (Cat. No. 555493/561889) antibodies. The erythrocytes were lysed with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899). The two-color flow cytometric dot plot shows the correlated expression of CD45RA versus CD45RO for gated events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

**Preparation and Storage**

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The antibody was conjugated with BD Horizon™ BV510 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV510 were removed.

**Application Notes****Application**

Flow cytometry	Routinely Tested
----------------	------------------

**Suggested Companion Products**

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 ml	(none)
563025	BV510 Mouse IgG2b, κ Isotype Control	50 µg	27-35
555899	Lysing Buffer	100 ml	(none)
555493	PE Mouse Anti-Human CD45RO	100 tests	UCHL1
561889	PE Mouse Anti-Human CD45RO	25 tests	UCHL1

**BD Biosciences**

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



## Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100- $\mu$ l experimental sample (a test).
2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
3. An isotype control should be used at the same concentration as the antibody of interest.
4. Please refer to [www.bdbiosciences.com/pharming/en/protocols](http://www.bdbiosciences.com/pharming/en/protocols) for technical protocols.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
7. Brilliant Violet™ 510 is a trademark of Sirigen.

## References

Knapp W, Dörken B, Gilks WR, et al, ed. *Leucocyte Typing IV*. New York, NY: Oxford University Press; 1989:1-1182. (Clone-specific)

Koristka S, Cartellieri M, Theil A, et al. Retargeting of human regulatory T cells by single-chain bispecific antibodies. *J Immunol*. 2012; 188(3):1551-1558.

(Clone-specific: Flow cytometry, Fluorescence activated cell sorting)

Zola H, Swart B, Nicholson I, Voss E. *Leukocyte and Stromal Cell Molecules. The CD Markers*. Hoboken, New Jersey: John Wiley & Sons, Inc.; 2007:1-581.

(Biology)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

