


Table of Contents	Page
1. Description	1
2. Overview of GeneBLAzer [®] GPCR Cell-based Assays	2
3. Overview of Division-Arrested Cells	2
4. Materials Supplied	3
5. Materials Required	4
5.1 Optional Equipment and Materials	4
6. Thawing DA Cells	5
7. Media Requirements	5
8. Assay Procedure	6
8.1 Quick Assay Reference Guides.....	6
8.2 Detailed Assay Protocol	8
8.3 Substrate Preparation, Loading and Incubation	9
8.4 Detection	9
9. Data Analysis.....	10
9.1 Background Subtraction and Ratio Calculation	10
9.2 Visual Observation of Intracellular Beta-lactamase Activity Using LiveBLAzer [™] -FRET B/G Substrate (CCF4-AM).....	10
10. References	11
11. Purchaser Notification	12

1. Description

GeneBLAzer[®] ADORA2A CHO-K1 DA (Division-arrested) cells contain the human Adenosine A2a Receptor (ADORA2A) stably integrated into the CellSensor[®] CRE-bla CHO-K1 cell line. CellSensor[®] CRE-bla CHO-K1 contains a beta-lactamase reporter gene under control of a CRE response element stably integrated into CHO-K1 cells. ADORA2A CHO-K1 DA cells have been functionally validated for Z'-Factor and EC₅₀ concentrations of 5'-(N-Ethylcarboxamido) adenosine (NECA).

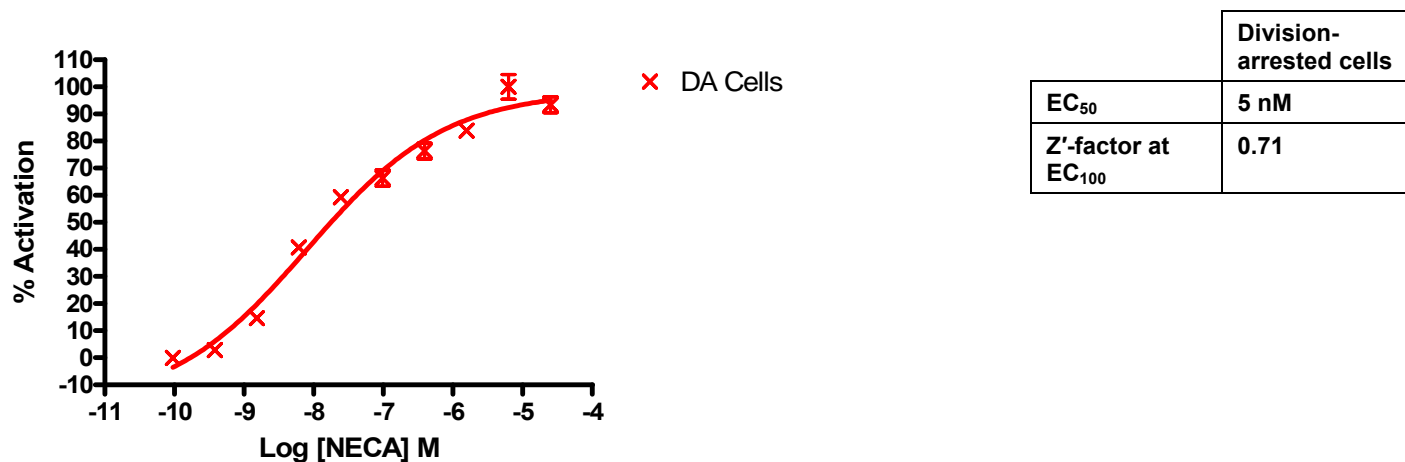


Dose response comparing ADORA2A CHO-K1 DA cells and ADORA2A-CRE-bla CHO-K1 cells to 5'-(N-Ethylcarboxamido) adenosine (NECA).

2. Overview of GeneBLAzer® GPCR Cell-based Assays

The GeneBLAzer® GPCR Cell-based Assay provides a highly accurate, sensitive, and easy-to-use method of monitoring cellular response to drug candidates or other stimuli (1). The core of the GeneBLAzer® assay technology is a beta-lactamase (*bla*) fluorescence resonance energy transfer (FRET) substrate that generates a ratiometric reporter response with minimal experimental noise. In addition to the dual-color (blue/green) readout of stimulated and unstimulated cells, this ratiometric method reduces the absolute and relative errors that can mask the underlying biological response of interest. Such errors include variations in cell number, transfection efficiency, substrate concentration, excitation path length, fluorescence detectors, and volume changes. The GeneBLAzer® GPCR assay technology has been proven in high-throughput screening (HTS) campaigns for a range of target classes, including G-protein coupled receptors (GPCRs) (2, 3), nuclear receptors (4-6) and kinase signaling pathways (7). The utility of division-arrested cells in HTS has also been demonstrated (8-11).

The GeneBLAzer® assay technology uses a mammalian-optimized *bla* reporter gene combined with a FRET-enabled substrate to provide reliable and sensitive detection in cells (1) (Figure 1). Cells are loaded with an engineered fluorescent substrate containing two fluorophores: coumarin and fluorescein. In the absence of *bla* expression, the substrate molecule remains intact. In this state, excitation of the coumarin results in fluorescence resonance energy transfer to the fluorescein moiety and emission of green fluorescent light. However, in the presence of *bla* expression, the substrate is cleaved separating the fluorophores and disrupting energy transfer. Excitation of the coumarin in the presence of *bla* enzyme activity results in a blue fluorescence signal.

Figure 1.

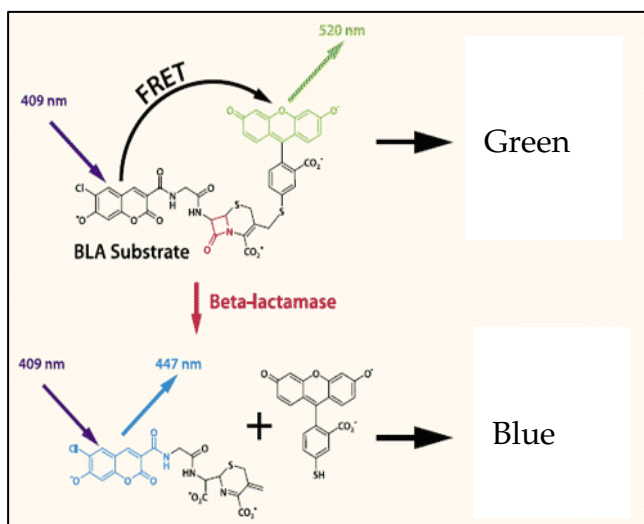

Fluorescent detection of beta-lactamase reporter gene response using a FRET-enabled substrate. After substrate loading, in the absence of beta-lactamase expression, cells appear green fluorescent. In the presence of beta-lactamase expression, the substrate is cleaved and cells appear blue fluorescent.

3. Overview of Division-Arrested Cells

Many of Invitrogen's cell lines are available in two forms: dividing or division arrested. Invitrogen's division-arrest technology allows the use of frozen cells, made from the exact same cell line sold in its dividing form, as ordinary, cost-effective assay reagents for screening. Division-arrested (DA) cells exhibit response profiles, similar to dividing cells thus ensuring you obtain the correct pharmacological profile.

DA cells may be plated and assayed within 24 hours of thawing. Cell numbers for DA cells increase only marginally after plating, thereby removing the variability caused by cell division during the course of an assay and providing more consistent results.

4. Materials Supplied

Product: (One of these)	Name	Size	Catalog #
	GeneBLAzer® ADORA2A CHO-K1 DA Assay Kit Each system contains sufficient division-arrested cells and substrate to assay one 384-well plate. Includes: <ul style="list-style-type: none">• GeneBLAzer® ADORA2A CHO-K1 DA cells (K1345A)• LiveBLAzer™-FRET B/G Loading Kit, 70 µg• Solution D, 1 ml• Protocol• Certificate of Analysis	1 plate	K1345
Shipping Condition:	Dry ice		
Storage Condition of Cells:	Liquid nitrogen. Immediately upon receipt, cells must be stored in liquid nitrogen or thawed for immediate use. Cells stored at -80°C can quickly lose viability.		
Cell Phenotype:	Epithelial		
Mycoplasma Testing:	Negative		
BioSafety Level:	1		

For Technical Support for this or other Drug Discovery Products, dial 760-603-7200, option 3, extension 40266

5. Materials Required

Use the table below to determine the additional media and reagents required for use with each kit:

Media/Reagents	Recommended Source	Part #	Required Separately?
			ADORA2A CHO-K1 DA Assay Kit (K1345)
LiveBLAzer [™] -FRET B/G Loading Kit: LiveBLAzer [™] -FRET B/G Substrate (CCF4-AM) DMSO for Solution A Solution B Solution C	Invitrogen	K1427 (70 µg) K1095 (200 µg) K1096 (1 mg) K1030 (5 mg)	No (included in kit)
Solution D	Invitrogen	K1156 (1 ml) K1157 (25 ml)	No (included in kit)
DMEM (high-glucose), with GlutaMAX [™]	Invitrogen	10569-010	Yes
DMSO	Fluka	41647	Yes
Fetal bovine serum (FBS), dialyzed, (DO NOT SUBSTITUTE!)	Invitrogen	26400-036	Yes
Non-essential amino acids (NEAA)	Invitrogen	11140-050	Yes
Penicillin/Streptomycin (antibiotics)	Invitrogen	15140-122	Yes
HEPES (1 M, pH 7.3)	Invitrogen	15630-080	Yes
5'-(N-Ethylcarboxamido) adenosine (NECA)	Sigma	E2387	Yes

The following table lists additional items required for use with all kits:

Consumables	Recommended Source	Part #
Black-wall, clear-bottom, 384-well assay plates (with low fluorescence background)	Corning	3712
Compressed air	Various	---
Equipment	Recommended Source	
Fluorescence plate reader with bottom-read capabilities	Various	
Filters if required for plate reader (see Section 8.4.1)	Chroma Technologies	

5.1 Optional Equipment and Materials

- Epifluorescence- or fluorescence-equipped microscope, with appropriate filters
- Microplate centrifuge

6. Thawing DA Cells

Note: Refer to **Section 7, Media Requirements** for specific media recipes.

Note: Once cells are thawed per the instructions below, cells must be counted and the density adjusted to the appropriate level as specified in **Section 8.2.2, Assay Procedure**, prior to analysis.

1. Rapidly thaw the vial of cells by placing at 37°C in a water bath with gentle agitation for 1–2 minutes. Do not submerge vial in water.
2. Decontaminate the vial by wiping with 70% ethanol before opening in a Class II biological safety cabinet.
3. Transfer the vial contents drop-wise into 10 ml of Assay Medium in a sterile 15-ml conical tube.
4. Centrifuge cells at 200 × *g* for 5 minutes.
5. Aspirate supernatant and resuspend the cell pellet in 1 ml fresh Assay Medium. (For the 1 × 384-well plate size of cells, dilute in 1 ml.)
6. Count the cells.
7. Adjust the cell density with Assay Medium to the appropriate cell density as specified in Section 8.2.2. Proceed to **Section 8, Assay Procedure**, for guidance on using cells in an assay.

7. Media Requirements

Note: Unless otherwise stated, have all media and solutions at least at room temperature (we recommend 37°C for optimal performance) before adding to cells.

Note: Make **NO MEDIA SUBSTITUTIONS**, as these cell lines have been specifically validated for optimal assay performance with these media.

Component	Plating Medium	Assay Medium
DMEM	90%	100%
Dialyzed FBS (Do not substitute!)	10%	—
NEAA	0.1 mM	—
HEPES (pH 7.3)	25 mM	—
Penicillin (antibiotic)	100 U/ml	—
Streptomycin (antibiotic)	100 µg/ml	—

8. Assay Procedure

The following instructions outline the recommended procedure for determining activity of compounds as modulators of ADORA2A using LiveBLAzer[™]-FRET B/G Substrate as the readout. If alternative substrates are used (*e.g.*, ToxBLAzer[™] DualScreen or LyticBLAzer[™] Loading kits), follow the loading protocol provided with the product.

8.1 Quick Assay Reference Guides

For a more detailed assay protocol, see **Section 8.2**.

Agonist Assay Quick Reference Guide

	Unstimulated Wells	Stimulated Wells	Cell-free Wells	Test Compound Wells
Step 1 Plate cells, incubate	32 μ l cells in Plating Medium (10,000 cells/well)	32 μ l cells in Plating Medium (10,000 cells/well)	32 μ l Plating Medium (no cells)	32 μ l cells in Plating Medium 10,000 cells/well)
Incubate cells for 16-20 hrs. at 37°C / 5%CO ₂				
Step 2 Media Change	Remove Growth Medium and Replenish with Assay Medium.			
Step 3 Add Agonist or Test Compounds	8 μ l Assay Medium with 0.5% DMSO	8 μ l 5X agonist in Assay Medium with 0.5% DMSO	8 μ l Assay Medium with 0.5% DMSO	8 μ l 5X Test Compounds in 0.5% DMSO
Step 4 Incubate cells	Incubate in a humidified 37°C/5% CO ₂ incubator for 5 hours			
Step 5 Prepare 6X Substrate Mix	6 μ l of 1 mM LiveBLAzer [™] -FRET B/G (CCF4-AM) Substrate + 60 μ l of solution B, mix. Add 904 μ l of Solution C, mix. Add 30 μ l of Solution D, mix.			
Step 6 Add Substrate Mixture	8 μ l per well			
Step 7 Incubate Substrate Mix. + cells	2 hours at room temperature in the dark			
Step 8 Detect activity	See Section 7.4			
Step 9 Analyze data	See Section 8			

Antagonist Assay Quick Reference Guide

	Unstimulated Wells	Stimulated Wells	Antagonist Control Wells	Cell-free Wells	Test Compound Wells
Step 1 Plate cells, incubate	32 μ l cells in Plating Medium (10,000 cells/well)	32 μ l cells in Plating Medium (10,000 cells/well)	32 μ l cells in Plating Medium (10,000 cells/well)	32 μ l Plating Medium (no cells)	32 μ l cells in Plating Medium (10,000 cells/well)
Incubate cells for 16-20 hrs. at 37°C / 5%CO ₂					
Step 2 Media Change	Remove Plating Medium and Replenish with Assay Medium.				
Step 3 Add Antagonist or Test Compounds, incubate	4 μ l Assay Medium with 0.5% DMSO	4 μ l Assay Medium with 0.5% DMSO	4 μ l 10X antagonist in Assay Medium with 0.5% DMSO	4 μ l Assay Medium with 0.5% DMSO	4 μ l 10X Test Compounds in Assay Medium with 0.5% DMSO
Incubate plate with Antagonist for 30 minutes before proceeding					
Step 4 Add Agonist	4 μ l Assay Medium with 0.5% DMSO	4 μ l 10X agonist in Assay Medium with 0.5% DMSO	4 μ l 10X agonist in Assay Medium with 0.5% DMSO	4 μ l 10X agonist in Assay Medium with 0.5% DMSO	4 μ l 10X agonist in Assay Medium with 0.5% DMSO
Step 5 Incubate cells	Incubate in a humidified 37°C/5% CO ₂ incubator for 5 hours				
Step 6 Prepare 6X Substrate Mix	6 μ l of 1 mM LiveBLAzer™-FRET B/G (CCF4-AM) Substrate + 60 μ l of solution B, mix. Add 904 μ l of Solution C, mix. Add 30 μ l of Solution D, mix.				
Step 7 Add Substrate Mixture	8 μ l per well				
Step 8 Incubate Mixture	2 hours at room temperature in the dark				
Step 9 Detect activity	See Section 7.4				
Step 10 Analyze data	See Section 8				

For Technical Support for this or other Drug Discovery Products, dial 760-603-7200, option 3, extension 40266

8.2 Detailed Assay Protocol

Plate layouts and experimental outlines will vary; in screening mode, we recommend using at least three wells for each control: Unstimulated Control, Stimulated Control, and Cell-free Control.

Note: Some solvents may affect assay performance. Assess the effects of solvent before screening. The cell stimulation procedure described below is carried out in the presence of 0.5% DMSO to simulate the effect that a Test Compound's solvent might have on the assay. If you use other solvents and/or solvent concentrations, optimize the following assay conditions appropriately.

8.2.1 Precautions

- Work on a dust-free, clean surface. Always handle the 384-well, black-wall, clear-bottom assay plate by the sides; do not touch the clear bottom of the assay plate.
- If pipetting manually, you may need to centrifuge the plate briefly at room temperature (for 1 minute at $14 \times g$) after additions to ensure all assay components are on the bottom of the wells.

8.2.2 Plating Cells

1. Thaw DA cells and resuspend in Assay Medium to a density of 312,500 cells/ml.
2. Add 32 μ l per well of the Assay Medium to the Cell-free Control wells. Add 32 μ l per well of the cell suspension to the Test Compound wells, the Unstimulated Control wells, and Stimulated Control wells. Incubate cells at 37°C/ 5% CO₂ for 16-20 hours.
3. Remove Plating Medium and replenish with 32 μ l Assay Medium. Proceed to **Section 8.2.3** for an Agonist assay or **Section 8.2.4** for an Antagonist assay.

8.2.3 Agonist Assay Plate Setup

Note: This subsection provides directions for performing an Agonist assay. See **Section 8.2.4** for directions for performing an Antagonist assay.

1. Prepare a stock solution of 0.5% DMSO in Assay Medium.
2. Prepare a 5X stock of Test Compounds in Assay Medium with 0.5% DMSO.
3. Prepare a 5X stock of agonist in Assay Medium with 0.5% DMSO. We recommend running a dose response curve to determine the optimal concentration of the agonist solution.
4. Add 8 μ l of the stock solution of 0.5% DMSO in Assay Medium to the Unstimulated Control and Cell-free Control wells.
5. Add 8 μ l of the 5X stock solution of agonist to the Stimulated Control wells.
6. Add 8 μ l of the 5X stock of Test Compounds to the Test Compound wells.
7. Incubate the Agonist assay plate in a humidified 37°C/5% CO₂ incubator for 5 hours. Then proceed to **Section 8.3** for Substrate Loading and Incubation.

8.2.4 Antagonist Assay Plate Setup

Note: This subsection provides directions for performing an Antagonist assay. See **Section 8.2.3** for directions for performing an Agonist assay.

1. Prepare a stock solution of 0.5% DMSO in Assay Medium.
2. Prepare a 10X stock of Test Compounds in Assay Medium with 0.5% DMSO.
3. Prepare a 10X stock of agonist in Assay Medium with 0.5% DMSO. We recommend running a dose response curve to determine the optimal agonist concentration. For antagonist assays, we recommend stimulating cells with an agonist concentration in the EC₅₀-EC₈₀ range.
4. Prepare a 10X stock of antagonist in Assay Medium with 0.5% DMSO. We recommend running a dose response curve to determine the optimal inhibition concentration for the Antagonist solution.
5. Add 4 μ l of the 10X stock of Test Compounds to the Test Compound wells.
6. Add 4 μ l of the stock solution of 0.5% DMSO to the Stimulated Control wells, the Unstimulated Control wells, and the Cell-free Control wells.
7. Add 4 μ l of the 10X stock of antagonist in Assay Medium with 0.5% DMSO to the Antagonist Control wells.
8. If desired, incubate the Test Compounds with the cells humidified 37°C/5% CO₂ incubator before proceeding. Typically, a 30-minute incubation is sufficient.
9. Add 4 μ l of the 10X stock solution of agonist to the Test Compound wells, the Stimulated Control wells, and the Antagonist Control wells.

10. Add 4 µl of Assay Medium with 0.5% DMSO to the Unstimulated Control and Cell-free Control wells.
11. Incubate the Antagonist assay plate in a humidified 37°C/5% CO₂ incubator for 5 hours. Then proceed to **Section 8.3** for Substrate Loading and Incubation.

8.3 Substrate Preparation, Loading and Incubation

This protocol is designed for loading cells with LiveBLAzer[™]-FRET B/G Substrate Mixture (CCF4-AM) Substrate Mixture. If you use alternative substrates, follow the loading protocol provided with the substrate.

Prepare LiveBLAzer[™]-FRET B/G Substrate Mixture (CCF4-AM) Substrate Mixture and load cells in the absence of direct strong lighting. Turn off the light in the hood.

1. Prepare Solution A: 1 mM LiveBLAzer[™]-FRET B/G Substrate (CCF4-AM) Substrate Mixture in dry DMSO by adding 912 µl of DMSO per mg of dry substrate. Store the aliquots of the stock solution at -20°C until use. The molecular weight of the LiveBLAzer[™]-FRET B/G Substrate (CCF4-AM) is 1096 g/mol.
2. Prepare 6X Loading Solution:
 - a. Add 6 µl of Solution A to 60 µl of Solution B and vortex.
 - b. Add 904 µl of Solution C to the above solution and vortex.
 - c. Add 30 µl of Solution D to the above solution and vortex.
3. Remove assay plate from the humidified 37°C/5% CO₂ incubator.
Note: Handle the plate gently and do not touch the bottom.
4. Add 8 µl of the 6X Substrate Mixture to each well.
5. Cover the plate to protect it from light and evaporation.
6. Incubate at room temperature for 2 hours.

8.4 Detection

Make measurements at room temperature from the bottom of the wells, preferably in 384-well, black-wall, clear-bottom assay plates with low fluorescence background. Before reading the plate, remove dust from the bottom with compressed air.

8.4.1 Instrumentation, Filters, and Plates

- Fluorescence plate reader with bottom reading capabilities.
- Recommended filters for fluorescence plate reader:
Excitation filter: 409/20 nm
Emission filter: 460/40 nm
Emission filter: 530/30 nm

8.4.2 Reading an Assay Plate

1. Set the fluorescence plate reader to bottom-read mode with optimal gain and 5 reads.
2. Allow the lamp in the fluorescence plate reader to warm up for at least 10 minutes before making measurements.
3. Use the following filter selections:

	Scan 1	Scan 2
Purpose:	Measure fluorescence in the Blue channel	Measure FRET signal in the Green channel
Excitation filter:	409/20 nm	409/20 nm
Emission filter:	460/40 nm	530/30 nm

9. Data Analysis

9.1 Background Subtraction and Ratio Calculation

We recommend that you subtract the background for both emission channels (460 nm and 530 nm).

1. Use the assay plate layout to identify the location of the Cell-free Control wells. These Control wells are used for background subtraction.
2. Determine the average emission from the Cell-free Control wells at both 460 nm (Average Blue Background) and 530 nm (Average Green Background).
3. Subtract the Average Blue background from all of the Blue emission data.
4. Subtract the Average Green background from all of the Green emission data.
5. Calculate the Blue/Green Emission Ratio for each well, by dividing the background-subtracted Blue emission values by the background-subtracted Green emission values.

9.2 Visual Observation of Intracellular Beta-lactamase Activity Using LiveBLAzer[™]-FRET B/G Substrate (CCF4-AM)

Note: Microscopic visualization of cells will cause photobleaching. Always read the assay plate in the fluorescence plate reader before performing microscopic visualization.

An inverted microscope equipped for epifluorescence and with either a xenon or mercury excitation lamp may be used to view the LiveBLAzer[™]-FRET B/G Substrate (CCF4-AM) signal in cells. To visually inspect the cells, you will need a long-pass filter passing blue and green fluorescence light, so that your eye can visually identify whether the cells are fluorescing green or blue.

Recommended filter sets for observing beta-lactamase activity are described below and are available from Chroma Technologies (800-824-7662, www.chroma.com).

Chroma Set # 41031

Excitation filter:	HQ405/20x (405 ± 10)
Dichroic mirror:	425 DCXR
Emission filter:	HQ435LP (435 long-pass)

Filter sizes vary for specific microscopes and need to be specified when the filters are ordered. For epifluorescence microscopes, a long-pass dichroic mirror is needed to separate excitation and emission light and should be matched to the excitation filter (to maximally block the excitation light around 405 nm, yet allow good transmission of the emitted light).

10. References

1. Zlokarnik, G., *et al*, **Quantitation of Transcription and Clonal Selection of Single Living Cells with Beta-Lactamase as Reporter**, (1998) *Science*; **279**: p84-88.
2. Kunapuli P, Ransom R, Murphy K, Pettibone D, Kerby J, Grimwood S, Zuck P, Hodder P, Lacson R, Hoffman I, Inglese J, Strulovici B, **Development of an Intact Cell Reporter Gene Beta-lactamase Assay for G Protein-coupled Receptors**, (2003) *Analytical Biochem.*; **314**: p16-29.
3. Xing, H., Pollok, B., *et al*, **A Fluorescent Reporter Assay For The Detection of Ligands Acting Through G1 Protein-coupled Receptors**, (2000) *J. Receptor & Signal Transduction Research*; **20**: p189-210.
4. Qureshi, S., *et al*, **A One-Arm Homologous Recombination Approach for Developing Nuclear Receptor Assays in Somatic Cells**, (2003) *Assay and Drug Dev. Tech.*; **1**: p755-766.
5. Peekhaus, N. *et al*, **A Beta-Lactamase-Dependent Gal4-Estrogen Receptor Transactivation Assay for the Ultra-High Throughput Screening of Estrogen Receptor Agonists in a 3,456-Well Format**, (2003) *Assay and Drug Dev Tech*, **1**: p789-800.
6. Chin, J., *et al*, **Miniaturization of Cell-Based, Beta-Lactamase-Dependent FRET Assays to Ultra-High Throughput Formats to Identify Agonists of Human Liver X Receptors**, (2003) *Assay and Drug Dev. Tech.*; **1**: p777-787.
7. Whitney M, Rockenstein E, Cantin G, Knapp T, Zlokarnik G, Sanders P, Durick K, Craig FF, Negulescu PA., **A Genome-wide Functional Assay of Signal Transduction in Living Mammalian Cells**, (1998) *Nat. Biotechnol.*; **16**: p1329-1333.
8. Fursov, *et al*, (2005) **Improving Consistency of cell-based assays by using division-arrested cells**. *Assay & Drug Development Technologies*, **Vol 3**, pp 7-15
9. Kunapuli, *et al*, (2005) **Application of division arrest technology to cell-based HTS: comparison with frozen and fresh cells**. *Assay & Drug Development Technologies*, **Vol 3**, pp17-26.
10. Digan, *et al*, (2005) **Evaluation of division-arrested cells for cell-based high-throughput screening and profiling**. *J. Biomolecular Screening*. **Vol. 10**, pp 615 – 623.
11. Vasudevan, *et al*, (2005) **Improving high-content-screening assay performance by using division-arrested cells**. *Assay & Drug Development Technologies*, **Vol 3**, p.p.515.

11. Purchaser Notification

Limited Use Label License No. 51: Blastidicin and the Blastidicin Selection Marker

Blasticidin and the blasticidin resistance gene (bsd) are the subject of U.S. Patent No. 5,527,701 sold under patent license for research purposes only. For information on purchasing a license to this product for purposes other than research, contact Licensing Department, Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, California 92008. Phone (760) 603-7200. Fax (760) 602-6500.

Limited Use Label License No. 108: Lentiviral Technology

The Lentiviral Technology (based upon the lentikat™ system) is licensed from Cell Genesys, Inc., under U.S. Patent Nos. 5,834,256; 5,858,740; 5,994,136; 6,013,516; 6,051,427; 6,165,782 and 6,218,187 and corresponding patents and applications in other countries for internal research purposes only. Use of this technology for gene therapy applications or bioprocessing other than for non-human research use requires a license from Cell Genesys (Cell Genesys, Inc. 342 Lakeside Drive, Foster City, California 94404). The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity), including non-gene therapy research and target validation applications in laboratory animals.

Limited Use Label License No. 109: Retroviral Helper Lines

Retroviral helper cell lines are licensed from Wisconsin Alumni Research Foundation, under U.S. Patents and corresponding patents and applications in other countries for internal research purposes only. Use of these cell lines for Commercial Purposes requires a license from Life Technologies.

Limited Use Label License No. 150: GeneBLAzer® Technology

This product and/or its use is the subject of one or more of U.S. Patent Nos. 5,741,657, 5,955,604, 6,291,162, and 6,472,205 and foreign equivalents licensed to Life Technologies Corporation. The right to use this product for internal research specifically excludes the right to use this product to identify, discover, and profile compounds that act as a flavor, fragrance or taste-enhancers and modify a target identified in taste, olfaction, or pheromone detection, which compound does not require FDA approval of an NDA for claims of safety and efficacy. The right to use methods claimed in the foregoing patents with this product for research purposes can only be acquired by the use of GeneBLAzer® substrates purchased from Life Technologies Corporation or its authorized distributors.

Limited Use Label License No. 315: Division Arrested Cell Products

The purchaser of this product may not remove or obtain components from cells sold under this label license or use these cells in a manner inconsistent with protocols provided by Life Technologies Corporation for use of division arrested cells. The purchaser of this product may not attempt to expand division arrested cells and may only use the purchased amount of this product.

Limited Use Label License No. 358: Research Use Only

The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed expressly, by implication, or by estoppel. This product is for internal research purposes only and is not for use in commercial applications of any kind, including, without limitation, quality control and commercial services such as reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact outlicensing@lifetech.com or Out Licensing, Life Technologies, 5791 Van Allen Way, Carlsbad, California 92008.

Use of Genetically Modified Organisms (GMO)

Information for European Customers The ADORA2A CHO-K1 DA cell lines are genetically modified to contain pLenti-bsd/CRE-bla (YFP) and pcDNA6.2 V5/Hygro ADORA2A As a condition of sale, use of this product must be in accordance with all applicable local legislation and guidelines including EC Directive 90/219/EEC on the contained use of genetically modified organisms.

© 2008, 2010 Invitrogen Corporation. All rights reserved. Reproduction forbidden without permission