

# MOUSE anti-HUMAN CD278 (ICOS)

Publication No. MAN0006808 Rev. 1.01

Store at 2° to 8°C

Catalog No.	Form	Amount	Excitation	Peak Emission	
A14754	FITC	25 tests	488 nm	519 nm	
A14767	PE	25 tests	496 nm	578 nm	
A14773	APC	25 tests	650 nm	660 nm	
A14783	PE-Cy®7	25 tests	488 nm	767 nm	

### **Product Description**

The Mouse anti-Human CD278 (ICOS) Monoclonal Antibody (mAb) reacts with human CD278, a T cell specific activation molecule and a member of the CD28/CTLA-4 family. Human CD278 is composed of 27 kDa and 29 kDa chains. CD278 on activated T cells has potent co-stimulatory activity for T cell activation, and is required for humoral immune responses (particularly for memory B cell and plasma cell generation). The CD278 ligand, B7h/B7RP-1 is expressed on activated antigen presenting cells, and on a number of inflamed peripheral tissues.

## **Product Specifications**

Clonality: Monoclonal Host/Class: Mouse IgG

Reactivity: Human CD278 (ICOS)

Alternate Names: ICOS (Inducible co-stimulatory

molecule), H4, CRP-1, AILIM

Apparent MW: 55–60 kDa

Sequence Identity:HumanClone/PAD:ISA-3Isotype: $IgG_1$ 

**Lot:** See product label

## **Product Applications**

Applications reported for the Mouse anti-Human CD278 (ICOS) mAb include flow cytometry.

Because conditions may vary, it is recommended that each investigator determine the optimal amount of antibody to be used for each application.

#### Stability

When stored as instructed, expires one year from date of receipt unless otherwise indicated on product label.

### Storage and Handling

Store reagents at 2° to 8°C. If the reagent is being diluted, it is recommended that only the quantity to be used within one week be diluted. Cells should be analyzed within 18 hours of staining for best results.

Avoid light exposure with fluorochrome-conjugated antibodies. Use dim light during handling, incubation with cells, and prior to analysis.

#### Storage Buffer

An aqueous buffer with 0.09% sodium azide, which may contain carrier protein/stabilizer.

**Caution:** Sodium azide is an extremely toxic and dangerous compound particularly when combined with acids or metals. Properly dispose of solutions containing sodium azide.

#### **Product Documentation**

To obtain a Certificate of Analysis or Safety Data Sheets (SDSs), visit www.lifetechnologies.com/support.

#### Related Products

Product Name	Quantity	Catalog no.
AbC <sup>™</sup> Anti-Mouse Bead Kit	1 kit	A10344
AbC <sup>™</sup> anti-Rat/Hamster Bead Kit	1 kit	A10389
Protein A Agarose	5 mL	15918-014
Recombinant Protein G (rProtein G) Agarose	5 mL	15920-010

**Explanation of symbols** 

Symbol	Description	Symbol	Description		
REF	Catalogue Number	LOT	Batch code		
RUO	Research Use Only	IVD	In vitro diagnostic medical device		
X	Use by	1	Temperature limitation		
***	Manufacturer	EC REP	European Community authorised representative		
[-]	Without, does not contain	[+]	With, contains		
from Light:	Protect from light	Æ	Consult accompanying documents		
$\bigcap_i$	Directs the user to consult instructions for use (IFU), accompanying the product.				

### **Limited Product Warranty**

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at <a href="https://www.lifetechnologies.com/termsandconditions">www.lifetechnologies.com/termsandconditions</a>. If you have any questions, please contact Life Technologies at <a href="https://www.lifetechnologies.com/support">www.lifetechnologies.com/support</a>

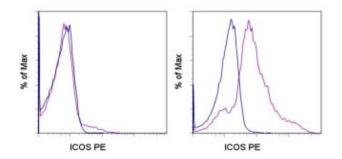
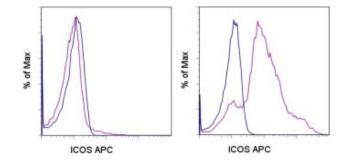


Figure 1 CD278 (ICOS) expression on human peripheral blood cells.

Staining of 3-day unstimulated (left) and 3-day anti-CD3/CD28-stimulated (right) human peripheral blood cells with a mouse IgG1 K-PE isotype control (blue histogram) or Mouse anti-Human CD278 (ICOS)-PE Monoclonal Antibody (Cat. no. A14767) (purple histogram). Cells in the lymphocytes gate were used for analysis.



**Figure 2** CD278 (ICOS) expression on human peripheral blood cells.

Staining of 3-day unstimulated (left) or 3-day anti-CD3/CD28-stimulated human peripheral blood cells with a mouse IgG1 K-APC isotype control (blue histogram) or Mouse anti-Human CD278 (ICOS)-APC Monoclonal Antibody (Cat. no. A14773) (purple histogram). Cells in the lymphocyte gate were used for analysis.

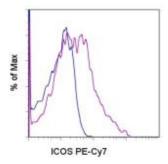


Figure 3 CD278 (ICOS) expression on human peripheral blood cells

Staining of normal human peripheral blood cells either unstimulated (blue histogram) or stimulated for 3 days with anti-CD3/CD28 (purple histogram) with Mouse anti-Human CD278 (ICOS)-PE-Cy®7 Monoclonal Antibody (Cat. no. A14783). Total viable cells were used for analysis.

**Note:** All flow cytometric data shown may not necessarily have been generated using the enclosed lot of reagent. For this reason, and due to differences in flow cytometers and cytometer settings, results may vary from those illustrated above. It is suggested that investigators titrate reagents to determine optimal conditions for use in their systems.

## Limited Use Label License: Research Use Only

The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed expressly, by implication, or by estoppel. This product is for internal research purposes only and is not for use in commercial applications of any kind, including, without limitation, quality control and commercial services such as reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact outlicensing@lifetech.com or Out Licensing, Life Technologies, 5791 Van Allen Way, Carlsbad, California 92008.

LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

©2012 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners. CV® is a trademark of

GE/Amersham Biosciences.

For support visit www.lifetechnologies.com/support or email techsupport@lifetech.com

