

# Bcl-xL Mouse Anti-Human Monoclonal Antibody (clone 7B2.5), Biotin Conjugate

Store at 2°C to 8°C

Pub. No. MAN0007963 Rev. 1.00

Catalog No.	Form	Amount	Excitation	Peak Emission
A15375	Biotin	0.1 mg	—	_

### **Product Description**

Apoptosis, or programmed cell death occurs in many cellular systems,<sup>1</sup> and plays a key role in tissue and organ development, and in cell turnover in adult tissues. Apoptosis can be induced by a variety of internal and external stimuli including growth factor deprivation, cytokine treatment, antigen-receptor engagement, cell-cell interactions, irradiation, and glucocorticoid treatment.<sup>2</sup> Bcl-2 and its homolog, Bcl-xL, protect cells from apoptosis,<sup>3,4</sup> while other homologs of Bcl-2 such as Bax, Bad, and Bak enhance apoptosis.<sup>5-8</sup> Bcl-xL blocks apoptosis induced by a variety of stimuli, and under certain conditions, offers greater protection against apoptosis than Bcl-2.9-13 In contrast, Bad and Bax inhibit the protective functions of Bcl-xL and Bcl-2, respectively. Although heterodimerization between Bcl xL/Bad and Bcl-2/Bax was originally thought to be essential for the differential anti-apoptotic activity of Bcl-xL and Bcl-2,<sup>5,14</sup> other results suggest that the formation of heterodimers may not be necessary for this death-repressing activity.<sup>15, 16</sup>

### **Product Specifications**

Clonality:	Μ
Host/Class:	М
Reactivity:	Н
Predicted Reactivity:	Н
Immunogen:	Re
Apparent MW:	29
Sequence Identity:	Н
Clone/PAD:	7E
Isotype:	Ig
Lot:	Se

Monoclonal Mouse IgG Human Bcl-xL Human, Rodent17 Recombinant Bcl-xS 29 kDa Human 7B2.5 IgG<sub>3</sub> See product label

### **Product Applications**

Applications reported for the Bcl-xL Mouse Anti-Human mAb include immunoprecipitation,<sup>15, 16</sup> immunohistochemistry,<sup>17</sup> flow cytometry,<sup>17</sup> and western blotting.<sup>18</sup>

Because conditions may vary, it is recommended that each investigator determine the optimal amount of antibody to be used for each application.

### Storage and Handling

Store reagents at 2°C to 8°C. If the reagent is being diluted, it is recommended that only the quantity to be used within one week be diluted. Cells should be analyzed within 18 hours of staining for best results.

Avoid light exposure with fluorochrome-conjugated antibodies. Use dim light during handling, incubation with cells, and prior to analysis.

### Stability

When stored as instructed, expires one year from date of receipt unless otherwise indicated on the Certificate of Analysis.

## Storage Buffer

Phosphate buffered saline (PBS) with 0.1% sodium azide.



**CAUTION!** Sodium azide is extremely toxic and may react with lead and copper plumbing to form highly explosive metal azides. Properly dispose of solutions containing sodium azide. Read the Safety Data Sheet (SDS) and follow the handling instructions. Wear appropriate protective eyewear, clothing, and gloves. SDSs are available at www.lifetechnologies.com/support.

### References


- 1. Cohen, J.J. 1991. Adv. Immunol. 50:55.
- 2. Cohen, J.J., and R.C. Duke. 1992. Annu. Rev. Immunol. 10:267.
- 3. Nunez, G., and M.F. Clarke. 1994. Trends Cell Biol. 4:399.
- 4. Cory, S. 1995. Annu. Rev. Immunol. 13:513.
- Oltvai., Z.N., C.L. Millman, and S.J. Korsmeyer. 1993. Cell 74:609.
- 6. Farrow, S.N., et al. 1995. Nature 374:731.
- 7. Chittenden, T., et al. 1995. Nature 374:733.
- 8. Kiefer, M.C., et al. 1995. Nature 374:736.
- Boise, L.H., M. Gonzalez-Garcia, C.E. Postema, L. Ding, T. Lindstein, L.A. Turka, X. Mao, G. Nunez, and C.B. Thompson 1993. *Cell* 74:597.
- 10. Gottschalk, A.R., L.H. Boise, C.B. Thompson, and J. Quintans. 1994. *Proc. Natl. Acad. Sci. USA* 91:7350.
- Gonzalez-Garcia, M., I. Garcia, L. Ding, S. O'Shea, L.H. Boise, C.B. Thompson, and G. Nunez. 1995. *Proc. Natl. Acad. Sci.* USA 92:4304
- 12. Dole, M.G., R. Jasty, M.J. Cooper, C.B. Thompson, G. Nunez, and V.P. Castle. 1995. *Can. Res.* 55:2576.
- 13. Shimizu, S., et al. 1995. Nature 374:811.
- 14. Yin, X.M., Z.N. Oltvai, and S.J. Korsmeyer. 1994. *Nature* 369:321.
- 15. Cheng, E.H.-Y., et al. 1996. Nature 379:554.
- 16. Gottschalk, A.R., L.H. Boise, Z.N. Oltvai, M.A. Accavitti, S.J. Korsmeyer, J. Quintans, and C.B. Thompson. 1996. *Cell Death and Differentiation* 3:113.
- 17. Boise, L.H. 1997. Personal communication.
- 18. L. Haughan. 1997. Personal communication.

#### **Product Documentation**

To obtain a Certificate of Analysis or Safety Data Sheets (SDSs), visit www.lifetechnologies.com/support.

#### Limited Product Warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at **www.lifetechnologies.com/termsandconditions**. If you have any questions, please contact Life Technologies at **www.lifetechnologies.com/support**.

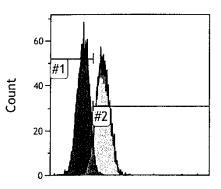

**Figure 1** Immunofluorescent staining of FL5 cells stained with Bcl-xL Mouse Anti-Human Monoclonal Antibody Biotin Conjugate (Cat. no. A15375) and FITC-conjugated Streptavidin.

Murine FL5 cells (#1) and FL5 cells transfected with Bcl xL expression plasmid (#2) were fixed with buffered paraformaldehyde and then permeabilized with saponin. The cells were incubated with Mouse Anti-Human Bcl-xL Monoclonal Antibody Biotin Conjugate (1 µg/10<sup>6</sup> cells) followed by FITC-conjugated Streptavidin.

### Related Products

Product Name	Quantity	Cat. No.			
AbC <sup>™</sup> Anti-Mouse Bead Kit	1 kit	A10344			
AbC <sup>™</sup> Anti-Rat/Hamster Bead Kit	1 kit	A10389			
FIX & PERM <sup>®</sup> Reagents (200 tests)	1 kit	GAS004			
LIVE/DEAD <sup>®</sup> Fixable Dead Cell Stain Kits					
Blue (UV excitation)		L23105			
Violet (405 nm excitation)	1 kit (200 assays)	L34955			
Yellow (405 nm excitation)		L34959			
Red (488 nm excitation)	(200 usbu ys)	L23102			
Near-IR (633/635 nm excitation)		L10119			

#### Explanation of symbols

Symbol	Description	Symbol	Description	Symbol	Description
***	Manufacturer	REF	Catalog number	LOT	Batch code
$\square$	Use by	X	Temperature limitation		
i	Consult instructions for use	$\triangle$	Caution, consult accompanying documents		

DISCLAIMER: LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

©2013 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation and/or its affiliates or their respective owners.

For support visit www.lifetechnologies.com/support or email techsupport@lifetech.com

