

Qdot®CD4 Rat Anti-Mouse Antibody Conjugates

Product	Form	Volume	Tests	Peak		Recommended Filter (nm)
				Excitation (nm)*	Emission (nm)	
Q10092	Qdot® 605	100 µL	100 minimum	405 (488)	605	605/20
Q22163	Qdot® 655	25 µL	25 minimum	405 (488)	655	660/20
Q22164	Qdot® 705	25 µL	25 minimum	405 (488)	705	695/40
Q22165	Qdot® 800	25 µL	25 minimum	405 (488)	800	780/60

*Qdot® nanocrystals excite optimally in the UV to 405 nm range, but can also be excited with wavelengths shorter than their emission maximum, such as with a 488-nm laser.

Product Description

Rat monoclonal antibody to the mouse CD4 antigen.

Product Specifications

Clone: RM4-5
Clonality: Monoclonal
Isotype: Rat IgG2a
Lot No.: See product label
Buffer: 50 mM borate, 1 M betaine, pH 8.3
Preservative: 0.05% sodium azide. Sodium azide is an extremely toxic and dangerous compound particularly when combined with acids or metals. Properly dispose of solutions containing sodium azide.

Storage and Handling

Store reagents at 2°C –8°C. **Do not freeze.** Because Qdot® nanocrystals are conjugated to biological materials, some loss of activity may be observed with prolonged storage.

Qdot® Antibody (Ab) conjugates are photostable, and do not need to be protected from light. However, if using Qdot® Ab conjugates in combination with conventional fluorochrome conjugated antibodies, minimize light exposure during handling, incubation with cells, and prior to analysis. We recommend analysis of cells within 18 hours of staining. If dilute reagent is used, dilute only the quantity of reagent to be used within one day.

The Qdot® Ab conjugates contain cadmium and selenium in an inorganic crystalline form. Dispose of the material in compliance with all applicable local, state, and federal regulations for disposal of these classes of material. For more information on the composition of these materials, consult the Safety Data Sheets (SDSs).

Stability

When stored as instructed, expires six months from date of receipt unless otherwise indicated on product label.

Qdot® Primary Antibody Conjugates

Qdot® Ab conjugates possess a bright fluorescence emission that makes them well suited for the detection of low-abundance extracellular proteins. Approximately the same size as R-phycoerythrin (R-PE) and compatible with existing organic fluorophore conjugates, Qdot® Ab conjugates can be excited with any wavelength below their emission maximum, but are best excited by UV or violet light. The narrow, symmetric emission profiles of Qdot® Ab conjugates allow for minimal compensation when using a single excitation source, and the very long stoke shifts enable better, more efficient multicolor assays using the 405 nm violet laser. Available in multiple colors for use flow cytometry, these advantages make Qdot® Ab conjugates powerful tools for antibody labeling and staining.^{1,2}

Product Characterization

Antigen Specificity: The RM4-5 monoclonal antibody (mAb) reacts with CD4 antigen which is expressed on the majority of thymocytes and on the MHC class II restricted subset of mature T cells, including Th cells.^{3,4} Mouse CD4 antigen has also been reported to be present on multipotential hematopoietic stem cells, bone marrow myeloid precursors, and intrathymic precursors.^{4,5} As a co-receptor in the TCR complex, CD4 is involved in T-cell activation through interaction with MHC class II on antigen presenting cells and in signal transduction via protein tyrosine kinase lck.³ The RM4-5 mAb can be used in immunostaining for flow cytometry and for immunohistochemistry on acetone-fixed, frozen tissue sections.

Leukocyte Workshop Status: NA

Product Use

Staining: Stain cells in any standard staining buffer, such as phosphate-buffered saline (PBS) with 1% bovine serum albumin (BSA). We recommend titrating the antibody conjugate to determine the optimal conditions for use in your specific system. Qdot® Ab conjugates may be mixed with other antibodies, but need to use diluted conjugates on day of dilution. Qdot® Ab conjugates can be used for surface staining applications with most conventional sample preparation reagents, such as Cal-Lyse™ Lysing Solution and FIX & PERM® reagents, with minimal effect on fluorescence. We have observed some batches of BD FACS® Lysing Solution to interfere with Qdot® Ab conjugate fluorescence.

Instrument setup: Qdot® Ab conjugates are excited optimally with UV or 405 nm light, although excitation can be obtained with any wavelength below the emission wavelength of a given nanocrystal. Make sure the cytometer has an appropriate emission filter for the Qdot® Ab conjugate being used. The table above has filter recommendations; alternate filters can be used as long as they capture the emission maximum, but filter width impacts spectral overlap corrections.

Note: Qdot® Ab conjugate can be used on cytometers that do not have UV or violet excitation sources as long as they have appropriate emission filters. Be sure to check for Qdot® Ab conjugate emission in any channel that can capture nanocrystal emission off of other lasers on the cytometer.

Product Quality Control

Each lot has been tested by flow cytometry using human peripheral blood leukocytes. See reverse for representative flow cytometry data.

For Research Use Only. Not for use in diagnostic procedures.

Life Technologies Corporation • 5791 Van Allen Way • Carlsbad • CA 92008 • Tel: 800.955.6288 • E-mail: techsupport@lifetech.com

References


1. Telford, W. G. 2004. Cytometry Part A 61A:9.
2. Perfetto, S. P., P. K. Chattopadhyah, and M. Roederer. 2004. Nature Reviews - Immunology 4: 648.
3. Bierer, B. E., B. P. Sleckman, S. E. Ratnofsky, and S. J. Burakoff. 1989. Ann. Rev. Immunol. 7: 579.
4. Fredrickson, G. G., and R. S. Basch. 1989. J. Exp. Med. 169: 1473.
5. Wu, L., R. Scollay, M. Egerton, M. Pearse, G. J. Spangrude, and K. Shortman. 1991. Nature 349: 71.

Related Products


Catalog no.	Product Name	Unit Size
GAS-010	Cal-Lyse™ Whole Blood Lysing Solution	25 mL
GAS-010S-100	Cal-Lyse™ Whole Blood Lysing Solution	100 mL
HYL-250	High-Yield Lyse Fixative	500 mL
GAS001S-5	FIX & PERM® Reagent A (Individual)	5 mL
GAS001S-100	FIX & PERM® Reagent A (Bulk)	100 mL
GAS002S-5	FIX & PERM® Reagent B (Individual)	5 mL
GAS002S-100	FIX & PERM® Reagent B (Bulk)	100 mL

Limited Product Warranty


Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at www.lifetechnologies.com/termsandconditions. If you have any questions, please contact Life Technologies at www.lifetechnologies.com/support.

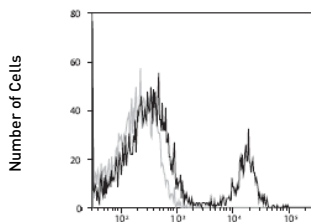


Qdot® 605 anti-mouse CD4 conjugate
605/20 band pass filter

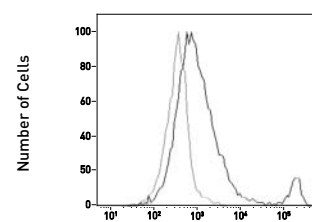

Qdot® 655 anti-mouse CD4 conjugate
660/20 band pass filter

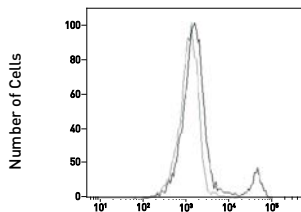

Qdot® 705 anti-mouse CD4 conjugate
695/40 band pass filter

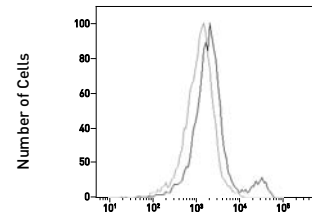

Qdot® 800 anti-mouse CD4 conjugate
780/60 band pass filter

Histograms of gated mouse splenocytes labeled with CD4 Anti-Mouse splenocytes mAb (clone S5.5) Qdot® Conjugates. Samples were acquired and analyzed using 405 nm excitation specified band pass emission filter on a BD LSR II flow cytometer (BD Biosciences, San Jose, CA) for the Qdot® 605 conjugate, or on an Attune® Acoustic Focusing Cytometer (Life Technologies, Carlsbad, CA) for the Qdot® 655, 705 and 800 conjugates. The black line histograms represent cells stained with anti-human CD4 antibody conjugate, and the gray line represents unstained cells.

Note: Flow cytometric data shown may not necessarily have been generated using the enclosed lot of reagent. For this reason, and due to differences in flow cytometers and cytometer settings, results may vary from those illustrated above. We recommend titrating reagents to determine optimal conditions for use in your systems.

Optimize your research with the Attune® Acoustic Focusing cytometer. Visit www.lifetechnologies.com/attune



DISCLAIMER: LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

©2012 Life Technologies Corporation. All rights reserved. This information is subject to change without notice.

The trademarks mentioned herein are the property of Life Technologies Corporation and/or its affiliates or their respective owners.

Fix & Perm is a registered trademark of An Der Grub Bio Research GMBH. BD FACS is a registered trademark of Becton, Dickinson and Company.



For support visit lifetechnologies.com/support or email techsupport@lifetech.com

lifetechnologies.com

11 January 2013